

SHOPRIDER®



Mikra

Model – TE-7A

USER MANUAL

Please ensure that this manual is read and understood before using the scooter.

Contents	Page
Introduction	3
Intended Use	3
Technical Specification	3
Feature Guide	4
Safety Advice	5
Adjustments	
• Tiller angle	6
• Seat swivel	6
• Seat removal	7
• Seat height	7
• Arm width adjustment	7
Tiller Controls	8
Brake Release / Freewheel Operation	9
Starting to Drive	
• Getting Started	10
• Moving Off	10
• Steering	10
• Reversing	10
• Ramps / Slopes	11
• Surfaces	11
• Transferring	11
Storage	11
Batteries and Battery Charger	11 - 12
Dismantling	13
Tiller Lock Lever	14
Maintenance	14
Troubleshooting	
• Fuses	15
• Check List	15
Further Information	16
Warranty	

Introduction

Congratulations on your purchase of the SHOPRIDER *Mikra* scooter. This scooter incorporates the latest innovative designs. The SHOPRIDER *Mikra* can be used in confined spaces both indoors and out. A direct drive rear transaxle and electromagnetic brake give a safe and efficient ride.

Intended Use

The SHOPRIDER *Mikra* is intended for indoor or outdoor use in pedestrian areas by a single individual with limited mobility up to a maximum weight of 113kg (18st). It is intended to be used by those suffering from some disability or impairment which prevents a normal level of personal mobility.

The acute manoeuvrability prohibits the use by persons who are severely physically impaired or those who have balance difficulties. Your assessor / healthcare professional will determine the suitability of the scooter for each individuals' intended use.

Note: Persons with lower body amputations should not use this device unless evaluated with artificial limbs fitted, as stability may be compromised due to an irregular centre of gravity.

Technical Specification

Class 1 Medical Device

Class 2 Indoor/Outdoor Vehicle

<i>Overall Size</i>		<i>Seat Dimensions</i>	
<i>Length</i>	1070mm / 42"	<i>Width</i>	470mm / 18"
<i>Width</i>	530mm / 21"	<i>Depth</i>	410mm / 16"
<i>Height (Max to top of seat)</i>	860mm / 34"	<i>Seat to floor pan</i>	365mm / 14" 390mm / 15"
<i>Battery</i>	2 x 12Ah	<i>Backrest height (Incl headrest)</i>	360mm / 14"
<i>Height (folded) No Seat</i>	450mm / 18"	<i>Range</i>	14km / 9m*
<i>Maximum user weight</i>	113kg / 18st	<i>Seat to floor</i>	560mm / 22" 610mm / 24"
<i>Safe working slope</i>	6°/ 10%	<i>Battery charger (off board)</i>	2 Amp
<i>Weight with batteries</i>	48kg / 106lbs	<i>Turning Radius</i>	1150mm / 45"
<i>Max. Speed (up to)</i>	6kmph / 4mph	<i>Ground clearance</i>	70mm / 3"

* Passenger weight in excess of 100kg, rough ground conditions, low temperatures and battery condition can affect maximum range.

Front Chassis – 15kg
Seat – 10kg

Rear Chassis – 10kg
Battery Pack – 6.5kg (each)



1. Tiller control head
2. Detachable swivel padded seat with fold down back
3. Detachable basket (Capacity 5kg)
4. Right side battery box
5. Left side battery box
6. Anti-tip wheels
7. Puncture proof non-marking wheels
8. Tiller angle adjustment knob



DO:

- Read this manual carefully before using your scooter.
- Turn the power off when transferring to and from the scooter.
- Select a slow speed when going down gradients (max 6°) or over uneven ground.
- Keep your feet on the scooter whilst driving.
- Ensure the seat is in a fixed position facing forwards before driving.
- Ensure the tiller is secure.
- Ensure the batteries are fully charged before your journey.
- Avoid rough or soft terrain and long grass wherever possible.
- Be aware of other road users and pedestrians.
- Turn the power off when not in use.
- Follow the maintenance guide to ensure safe operation of the scooter.
- Always remove the keys if leaving the scooter unattended. Unauthorised persons/children may use the scooter resulting in injury.



DO NOT:

- Carry passengers.
- Exceed inclines over 6°.
- Traverse inclines.
- Attempt inclines if the surface is slippery.
- Turn sharply at full speed.
- Exit or enter the scooter unless the power is switched off.
- Switch off the power whilst driving – this will cause the scooter to stop abruptly (see emergency procedures – page 16).
- Use the scooter where you cannot safely or legally walk.
- Attempt to climb or descend kerbs except drop down kerbs.
- Use the scooter if under the influence of alcohol or any medication.
- Use a mobile phone whilst driving.
- Attempt to negotiate steps or escalators.
- Sit on the scooter if being transported in a vehicle.
- Attempt to tow any other devices.
- Exceed the stated weight limit.
- Connect any other devices to the batteries as this could cause permanent damage and void your warranty.
- Use any other batteries/charger other than those supplied by *SHOPRIDER*.
- Travel through water or over soft terrain (mud, sand etc.)
- Lift the scooter by the shroud assembly or any plastic part.
- Sharply shift your centre of gravity as this could make the scooter unstable.
- Freewheel down hills. This will damage the motor and brake and void the warranty,
- Tow the scooter with any other vehicle in the event of a breakdown.

Every care has been taken to ensure that your scooter will give you safe, reliable service provided you adhere to the advice given above.

Adjustments


Tiller Angle Adjustment

At the bottom of the tiller stem is the angle adjustment knob. Rotate the hand knob anticlockwise to loosen and adjust to the required angle. Re-tighten the hand knob fully before use. ***DO NOT fully undo the hand knob as the internal components may become detached from the assembly.***



Seat Swivel

The seat has a lockable swivel base to allow easier transfers. Simply lift the end of the lever and rotate the seat at the same time. Releasing the lever will allow the seat to lock in position. The seat locks in 45-degree intervals. Always ensure the seat is in a locked position before mounting or dis-mounting the scooter. The seat should always be locked facing forward whilst driving.

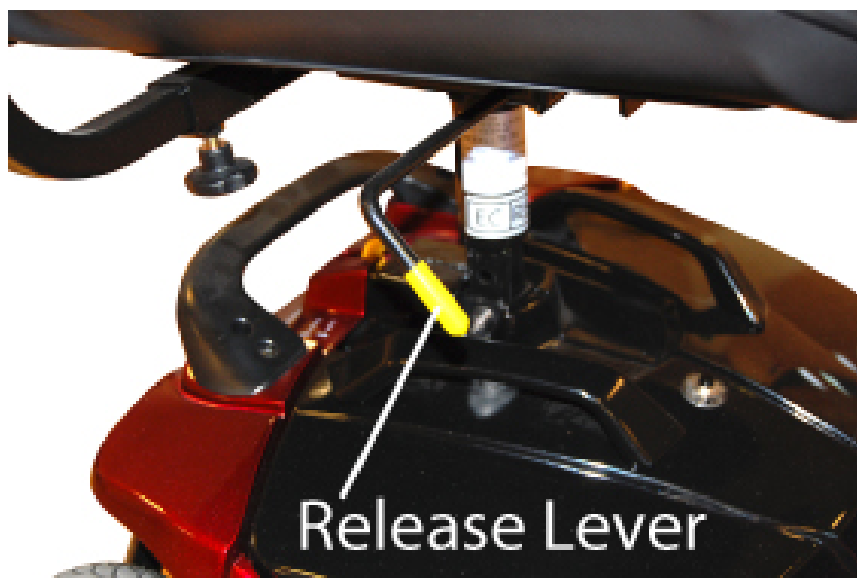

Seat Removal

To remove the seat, fold down the back of the seat and simply lift the seat away from the chassis while lifting the seat swivel lever.


Do not hang baggage or other objects on the back of the seat as this may affect stability.

Battery Box Removal

The battery boxes are simply lifted out of position. The LEFT battery box MUST be removed before the right hand side. The procedure is reversed for refitting the battery packs. (*Left and right sides are taken as if sitting on the scooter facing forward*).

Seat Height Adjustment

Remove the seat as previously described and remove the battery box as described above. Undo the stem-locking bolt and adjust the stem to the desired height. Re-fit the bolt and tighten. Re-fit the seat as above.

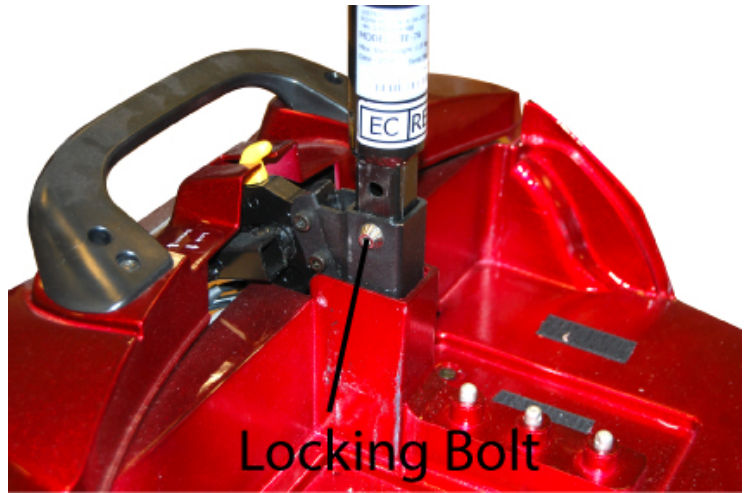


The higher the seat is adjusted; stability may be affected due to the higher centre of gravity.

Arm Adjustments

On each side of the seat towards the back there are two hand knobs, these are to adjust the width of the armrests or to remove them for transport.

Turn the hand knobs anti-clockwise to loosen. Once the knobs are loose you will be able to adjust the armrests in or out to the desired width. Re-tighten when a comfortable width has been achieved or remove the armrests completely for transport.

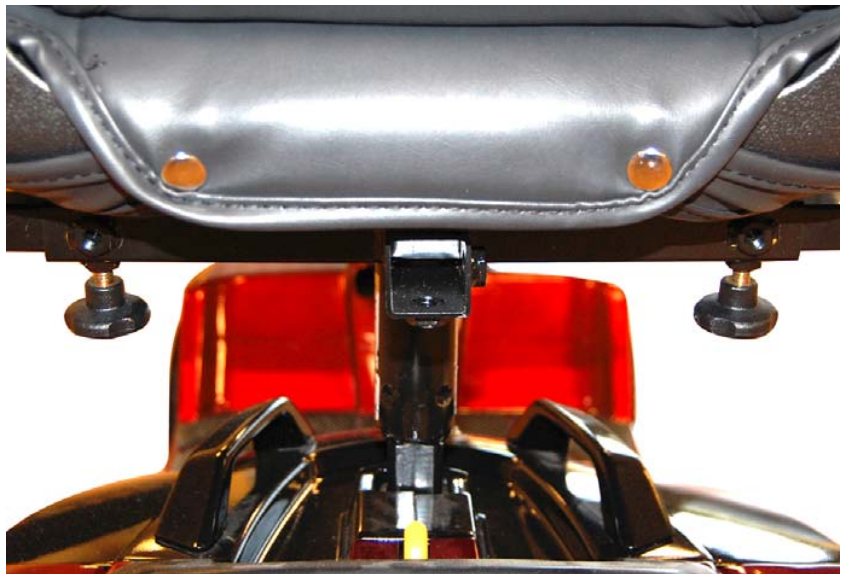

If removing the armrests for transport, ensure that the hand knobs are re-tightened to avoid them rattling loose and falling out.

Do not apply full body weight to the armrests for transferring. Lift the arm up out of the way to transfer to another seat.



Maximum Speed Knob

The speed control enables you to set the maximum speed of the scooter. Turn clockwise to increase the maximum speed and anti-clockwise to decrease. ***Always start at a slow speed until accustomed with the characteristics of your scooter. Do not put excessive pressure on the dial. Only turn as far as the stops.***

When driving in pedestrian areas such as shopping precincts and busy pavement areas always respect other pedestrians and turn the maximum speed down so that the scooter travels at the walking pace of others around you.

Horn Button

The yellow button operates the scooters horn. ***This should only be used to warn others of your presence and not as a rebuke of any form.***

Finger Lever (Wig-Wag)

The finger levers control the reverse and forward motion of the scooter – see 'Starting To Drive' for operation. Remember, more pressure = more speed. Release the lever and the scooter will automatically come to a safe stop as the electromagnetic brake is applied. ***Always release the lever if the scooter behaves erratically.***

On/Off Key

The key switch is located on the front of the tiller. Simply insert the key and turn clockwise to switch the power on. The display will light up to indicate the power is on. There may be a slight delay before the thumb levers operate as the controller checks the scooters circuits. ***Never leave the key in the switch unattended to prevent unauthorised use. Do not switch the power off to stop the scooter unless it is an emergency. The scooter will stop very abruptly if stopped in this way.***

Battery Level Indicator Gauge

This gauge gives you a guide to the amount of charge left in the scooters batteries. This will give the truest reading when the scooter is in motion. ***Always fully charge the batteries before driving your scooter and never let the batteries become completely drained, as this may cause permanent damage to the batteries.***

Brake Release / Freewheel Operation

The brake release lever is situated at the right rear of the scooter. This lever allows you to disengage the drive mechanism and push the scooter in 'freewheel' mode. To disengage the drive, pull the lever up as shown in the image below. The scooter can now be pushed.

Simply push the lever down toward the ground to re-engage the drive mechanism (the scooter can now be driven). The scooter will not operate when the lever is in the disengaged position. The scooter must be turned off and on by the key switch to reset the circuits before it will resume normal operation



Caution! Never leave the scooter unattended whilst in freewheel mode as it will roll away down an incline. 'Freewheel' mode should only be used in an emergency to move the scooter due to power failure. NEVER move the scooter in this mode whilst it is occupied. NEVER tow the scooter.

When in freewheel mode the scooter may be pushed on level ground at a walking pace. The scooter will 'auto-brake' if this speed is exceeded.

DO NOT step on the lever to re-engage the brake as this will damage the brake and micro-switch.

NEVER sit on the scooter and freewheel down a hill as this can cause permanent damage to the motor, brake and controller due to excessive regenerative braking. As above the scooter can only be moved at a walking pace when in freewheel mode.

Starting to Drive

Getting Started



Note: You must fully charge the batteries prior to first use (see 'Batteries and Charging')

Select an area with plenty of space.

Mount the scooter by lifting the arm and sliding onto the seat. **Ensure the seat is securely locked in the forward facing position before driving.** Adjust the tiller and seat if necessary (see 'Adjustments'). **DO NOT use the tiller to fully support you during transfers.**

The Mikra has a 'wig - wag' finger lever. The direction of the scooter is determined by which part of the lever is pushed / pulled as outlined above. Pulling the right side of the lever will cause the scooter to move forwards, the left side will select reverse.



Moving Off

Insert the key and switch on the power and initially set the speed control to a low setting. The further you pull the lever the faster the scooter will travel, but only to the maximum speed that has been selected. Release the lever and the scooter will come to a smooth controlled stop.



Caution! DO NOT switch off the power to stop the scooter as this will cause the transmission to lock and the scooter will come to a sharp stop, which could possibly cause injury. This procedure should only be used in an emergency.

Steering

Steering the scooter is light and easy. Ensure you give obstacles a wide clearance to avoid the rear wheels catching. When steering through limited spaces, turn the tiller in the direction that you wish to go and apply power gently. The scooter can turn very tightly with complete stability at slow speeds.



Caution! NEVER turn sharply at high speeds as it can make the scooter unstable. The maximum speed should be set at 50%.

Reversing



Before reversing ensure the area behind the scooter is free from any obstacles or hazards. Adjust the speed setting to low before you reverse. Steer wide of all corners and obstacles.

Caution! Do not reverse onto a main thoroughfare.

Ramps / Slopes

The stability of your scooter depends on several factors such as seat height, angle of gradient and the users' height and weight. When approaching an incline, do so directly and not at an angle and avoid making turns where possible.



Surfaces

Your scooter will perform well over many surfaces. However, avoid long grass, loose gravel, sand, mud and ice.

Transferring

Lift up the arm pad from the side you wish to leave or sit on the scooter. Use an aid for balance if required. DO NOT use the tiller or armrests to support the full weight during transfers.

Storage

Your scooter should be stored in a dry environment, preferably a garage or shed. Tarpaulins and covers may create condensation during certain weather conditions.

If your scooter is to be stored for prolonged periods, fully charge the batteries and remove them from the scooter.

Very cold weather can severely deplete battery condition (up to 40%). If left for any period during cold weather it is very important that the scooter is fully charged before use.

Batteries and Battery Chargers

Batteries

Charge the batteries completely (approx 8 hours) prior to first use.

The SHOPRIDER *Mikra* uses two 12Amp, 12 Volt sealed VRLA AGM batteries. The battery packs are located under the front of the seat. They can be easily removed by simply lifting by the handle on the top.

Low temperatures, rough terrain and the weight of the user can affect battery performance. Your battery gauge is a guide only to the level of charge remaining in your batteries and will give its best indication while the scooter is moving.

Note: Batteries recommended by SHOPRIDER MUST be used if they need to be replaced. An authorised SHOPRIDER dealer must carry out battery replacement.

Battery Charger

The battery charger is a 2-amp switch mode type off board. The battery charging ports are located on the left side of the tiller head and on the front of the battery pack. Before using your scooter for the first time be sure to fully charge the batteries, this may take up to 12 hours.

To charge the batteries connect the 3-pin male charger plug into the charging socket on the left side of the tiller. Plug the 3-pin mains plug into a wall socket and switch on. A constant yellow light will show when the battery charger is connected and charging is underway. The green LED will light when the charging cycle is complete. A red LED will light when there is a break in the connection. Always turn off mains power to the charger before disconnecting the plug from the battery box.



Caution! Do not charge the batteries for longer than 12 hours. Always charge the batteries in a well-ventilated indoor area. You must still maintain charging the batteries even if the scooter is not to be used for a prolonged period of time (2 to 3 weeks). Failure to do so could result in permanent damage and render the batteries unusable. Wherever possible do not interrupt the charge cycle until the batteries are fully charged. Continual interruption can permanently damage the batteries. If the charging lead should get damaged discontinue use immediately and seek a replacement from your dealer. **ONLY USE THE CHARGER SUPPLIED WITH YOUR SCOOTER.** Do not leave the charger plugged into to the charger port after the charge is complete. **ALWAYS** remove the keys from the scooter when charging.



Disposing and recycling of used batteries

Batteries must always be recycled correctly. Do not dispose of them with your domestic refuse.



For further information regarding your nearest recycling facility please contact your local authority.

Circuit Breaker



Battery Charging Socket

Dismantling

The scooter can be easily and quickly dismantled into five parts for transportation and storage. No tools are required.


Firstly remove the basket and seat as described earlier.

Remove the seat and battery packs (left side first).

Locate the release lever as shown below. The lever is sprung, push it to the rear of the scooter and lift the main body by the seat post as shown. The scooter will then separate.



To reassemble tilt the rear part of the scooter backwards and locate the rear frame tube into the frame locating lugs (shown right). The scooter should then click back together. Lift the release lever to ensure that the locking mechanism fully relocates and that there is no debris fouling the mechanism.

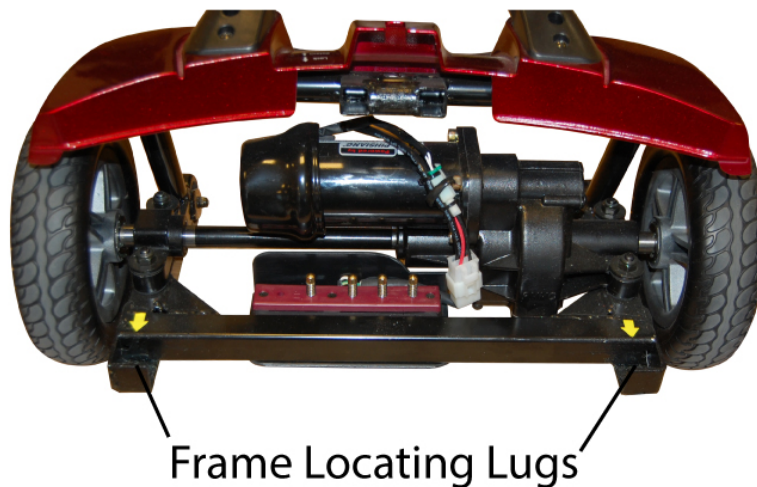


Tiller Lock Lever

Please refer to the image below. To lock the tiller direction the wheels should be in the forward direction. Lift the locking lever up and then slide it to the left as indicated by the arrows.

To unlock, simply slide the lever to the right.

CAUTION: *Never drive the scooter with the tiller in the locked position as this could lead to injury.*

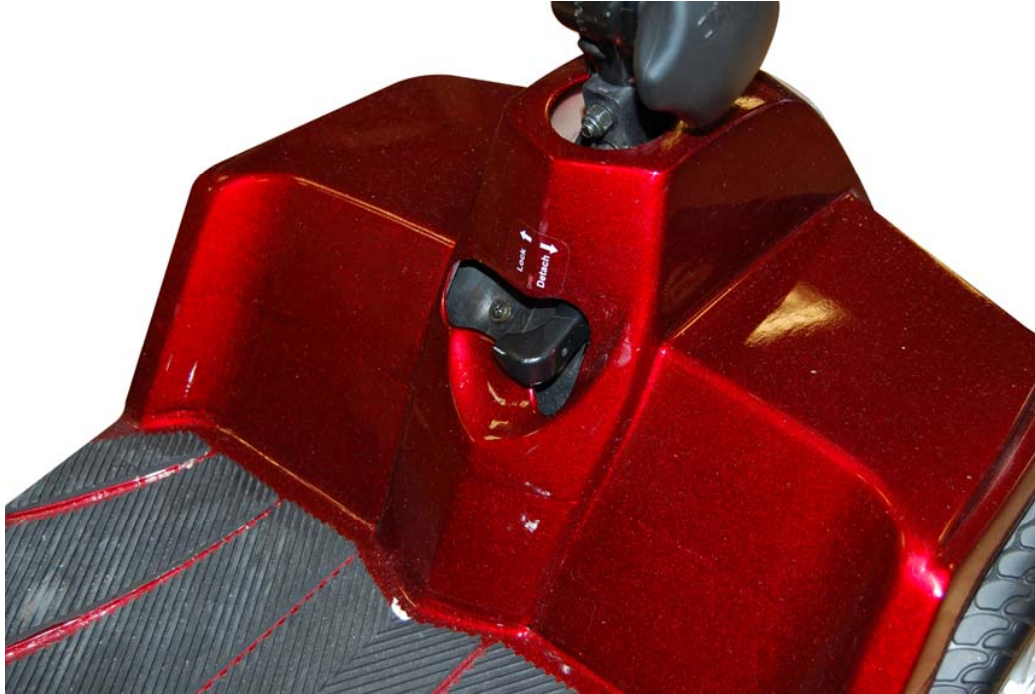

Maintenance

Product Life Expectancy

Your scooter has been designed to provide the user with a long, reliable life provided it is correctly maintained and regularly serviced by an authorised dealer. However, due to inevitable wear and tear and technological improvements it is recommended that the average useable life of this scooter is five years.

General Cleaning

The shrouds and seat of the scooter should be cleaned using a damp cloth and gentle detergent. Do not use abrasive cleaners or solvents, as this will damage the plastic components. Do not spray water (*hose or pressure washer*) onto the scooter as this may cause damage to the electronics.

Routine Maintenance

We suggest that your scooter should undergo a full service at least once a year by an authorised SHOPRIDER dealer.

Please check the following on a weekly basis:

- Front wheels spin freely
- Moving parts, pivots and frame fixings are not damaged or worn
- All hand knobs and levers are secure

Troubleshooting

Your scooter is fitted with the latest in electronic controls that are programmed to protect the electrical system from abnormal overloads.

Fuses

There are two fuses on the battery cables to give protection against possible overload on the wiring.

Check List

If your scooter fails to operate, please check the following:

- The unit is switched on.
- All plugs and connections are firmly fixed.
- Battery level shows full charge.
- Freewheeling device is in the drive position.
- Circuit breaker has not tripped.

If the freewheel device has been operated with the scooter switched on, the scooter will not operate until the key switch has been turned off and on again.

Should a fault occur, the fault code indicator on the tiller head will indicate a flash sequence which relates to the fault in the table below.



Safety Information Regarding Electromagnetic Interference (EMI)

Scooters and powerchairs are designed to operate under certain conditions. However, radio waves or electromagnetic fields may affect the operation of the scooter. The source of EMI could be radio or television transmitters or portable devices such as mobile phones, portable CB radios etc.

Should any unintentional movement occur turn the scooter off as soon as safely possible.

Solar Energy

During periods of sunny weather, it is possible that surfaces of your scooter may reach high temperatures (as with motor vehicles). Please take precautions if leaving the scooter in direct sunlight for any period of time. Seek shade or cover the seat and hand controls to keep any temperature increases to a minimum.

Health Factors

Suitability for this product has been assessed on the user's current state of health. It is essential that the user regularly monitor their ability to use the scooter safely.

General Use

Take heed of the information given in this manual regarding safety and the maintenance of your scooter. If adhered to, your scooter should give you years of reliable mobility.

Storage and General Care

We recommend that the scooter is always stored in a dry environment such as a garage or shed to protect the chassis and other components from corrosion. Do not leave the scooter outside during rainstorms. If the scooter is wet, wipe with a dry cloth to ensure there is no ingress of water into the internal components of the scooter.

During everyday use, do not use in rain showers or travel through puddles. Avoid use in areas where the scooter could become contaminated with sand or abrasive materials.

Emergency Procedures

Should your scooter behave erratically or become unresponsive whilst driving, release the wig-wag lever so it may return to its neutral position and switch off the power at the key / switch. These procedures will stop the scooter and activate the braking system, preventing any further uncontrolled manoeuvres.

Note: Switching the power off at the key / switch can cause the scooter to stop abruptly and should only be used to stop the scooter during an emergency as this can damage the drive unit.

Product Stability

Scooters and powered wheelchairs are designed to allow the user flexibility of movement in restricted areas. Driving at maximum speed and attempting sharp turns can create a situation where the product becomes unstable. Therefore it is essential that the user pay attention to the guidance set out in this manual.

For maximum stability the speed setting should be set up to 50% of the maximum available (halfway between the tortoise and the hare or just entering the red sector on the speed indicator depending on the model.)

The main factors affecting stability are:

- Maximum speed setting
- Seat height
- User movement
- User's manoeuvres
- Gradients or slopes

Please keep your speed to a level where you can safely stop without having to make any erratic turns.

Do not load bags etc. on the back of the seat as this will affect the centre of gravity and could cause instability.

Do not lean outside of the 'footprint' of the product, especially whilst moving.

Avoid making turns on gradients or slopes.

Driving in Pedestrian Areas

As mentioned previously, when using the scooter in pedestrian areas always lower your maximum speed to that of the pedestrians around you. Remember the stopping distance of your scooter and always allow enough space to stop safely to avoid collisions with pedestrians or other objects.

Warranty Terms and Conditions

This is to certify that your Shoprider® product is warranted for a period of 12 months from the date of original purchase. This warranty is not transferable.

Should a defect or operating fault arise within this period the Dealer from where the product was originally purchased should be notified immediately. Failure to register the fault may invalidate the warranty. The Dealer from where the product was purchased is responsible for any warranty repairs.

If a fault occurs within this period due to faulty workmanship or materials that the product will be repaired (or replaced at the Companys discretion) free of charge.

Any unauthorised modifications will invalidate the warranty.

Parts which are subjected to normal wear and tear, accidental damage, neglect or misuse are excluded from this warranty. Normal wear and tear includes tyres, chassis covers, bulbs, seat coverings etc.

Component parts replaced outside this period carry a six-month warranty. (Standard terms and conditions apply for this period).

Model: _____

Serial No: _____

Purchase Date: _____

Dealer Stamp

Distributed by:-
Roma Medical Aids Ltds
York Road, Bridgend Industrial Estate,
Bridgend
CF31 3TB
Tel: 01656 674488 Fax: 01656 674499
Email: sales@romamedical.co.uk
www.romamedical.co.uk

Shoprider reserve the right to change specifications without prior notice



Manufactured by:
Phisiang Machinery MFG. CO. LTD
No. 108 Hsin-He Rd, Hsin-Feng Hsiang, Hsinchu County, Taiwan

