

ICON 125

THE COMFORT WHEELCHAIR



REHASENSE

**NEW
PRODUCT**



ICON 125

GREAT SEATING AND POSITIONAL COMFORT



Comfort wheelchairs, or “multifunctional wheelchairs” as they are called, have to meet high requirements in order to meet the needs of users and caregivers. ICON 125 provide great seating and positional comfort as well as a lot of inbuilt adjustment capabilities to ensure a precise fit. Our modular ICON mobility system makes it possible to exchange and combine many components of our comfort wheelchairs ICON 120 (made of steel) and 125 (made of aluminium) with wheelchairs from our ICON 10-40 range.

FEATURES AT A GLANCE

- Hip widths can be adjusted with the arm rests and side panel pads
- The laminated seat and back cushions provide good postural control and pressure distribution over bony prominence with good lateral support to prevent the user from drifting sideways
- The back frame is angle, depth and height adjustable to fit all kinds of backs
- The head/neck support has 11,5 cm of medio-lateral slide adjustment as well as height, angle and depth adjustment
- The ICON Mobility System allows infinite adjustment of the seat depth

COLOURS

Silver

TECHNICAL DATA

SWL	Weight	Seat widths (SW)	Seat depths	Front seat height	Rear seat height	Tilt angle	Recline range
150 kg	SW 41 cm 31 kg SW 46 cm 33 kg SW 51 cm 36 kg	41, 46, 51 cm	40 – 53 cm (for seat width 41/46) 45 – 58 cm (for seat width 51)	44 – 47 cm	43 – 50.5 cm	up to 35 °	up to 30 °

WHEELS & BRAKES

- Can be equipped with 20”, 22” and 24” rear wheels for self-propelled users as well as with 16” rear wheels for attendants
- Front wheels are available in PU black with 6”, 7” and 8”
- Brake options include drum brakes and knee brakes with various grip options

FRAME & TRANSPORT

- Wheelbase can be extended at the front by 5 cm (SW: 41/46 cm) or 10 cm (SW: 51 cm)
- Equipped with front and rear mounting brackets for attaching the wheelchair as a seat in the car

BACK AND SEAT CUSHIONS

- Back: good lateral and lumbar support
- Seat: pressure relief under the ischial tuberosity, additional support on the lateral and medial sides of the thighs and in the rear area of the buttock rim
- Abduction pommel is also available

SPECIAL FEATURES




Arm rests: drop-in, height & depth adjustable
Leg rests with calf pads: swing away, detachable & liftable



Fully adjustable headrest, with folding side wings for more lateral support



Height adjustable push handles

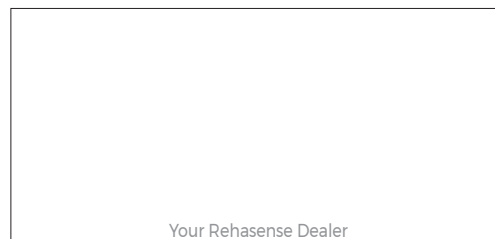

VI-2023

Rehasense's quality management complies with both the ISO 9001 and ISO 13485 standards, as confirmed by TÜV NORD. Rehasense's medical devices are CE-marked and meet the requirements of EU Regulation 2017/745 (MDR).

